

MEDIPHOT 937/943 X-RAY



COLENTA Medical X-Ray Film Processors: MP937 and MP943

For all standard X-Ray Films - Certified for **Mammography - Film** processing.

	DEV-Time	Cycle Time	Intake Speed
MP937/943:	18sec to	1'24''min to	106cm/min to (min.)
max. 37 (43) cm	99 sec	7'43''min	20cm/min (max.)

All **MEDIPHOT**-model processors have 9 programmable processing channels, an intermediate rinse bath rack system, film feed tray (optional with Cover lid) and supplied with 2 x 30 lt replenishment storage tanks.


VERSION I




Standard Features and accessories:

- compact / robust design
- open stand (optional covered)
- 2 rep. storage tanks 30L (for Dev and Fix)
- optical sensorbar for touchless measuring incoming film area

VERSION I Media Exit on Top Cover Lid

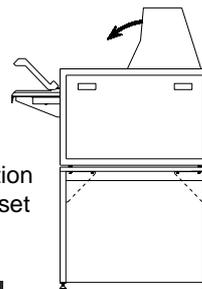

VERSION II Rear Media Exit in Receiving Basket

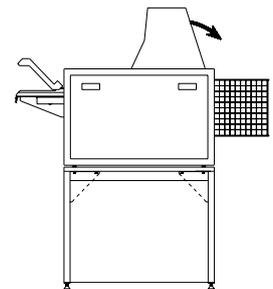


VERSION II

The best results in „Low Volume” Tank designs

All **Colenta** Processors are equipped with precision temperature control of all processing chemistry with automatic Developer Tank cooling to prevent overheating. Replenishment is preheated by injection into the inline chemical heating systems, which are controlled to pre-set programmed temperatures within +/- 0,2 °C.

TRANSPORT SYSTEM:

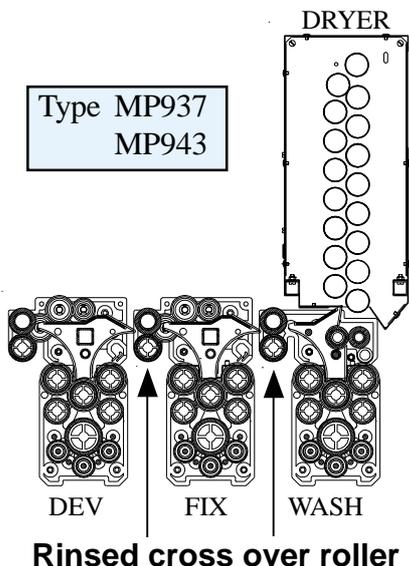

Optional accessories:

Rack Carrier Tray/Drip Tray for easy handling the transport racks.

Cover Lid for Film Feed Tray



Microprocessor controlled Displa



Closed stand for processor

Colenta®

Technical specification subject to change without notice.

see additional Data sheet: N1-15

COLENTA LABORTECHNIK
Ges.m.b.H. & Co. KG
Neunkirchner Strasse 117
A-2700 Wiener Neustadt
Tel.: +43-2622-28311-0
Fax: +43-2622-28311-7
E-mail: Colenta@nextra.at
<http://www.colenta.de>

ACI-COLENTA
INTERNATIONAL Div.
Forge Close, Little End Rd,
Eaton Socon, Cambs. UK
Tel. +44-1480-211118
Fax: +44-1480-215236
Fax: +44/870/132 1352
E-mail: aci@patflaherty.com